


Technical Specifications & Requirements for Software Installation and Use


Cardea SOLO 7-day Cardiac Monitor System

About Cardea SOLO® System

Cardea SOLO is the first complete system for timely, cost-effective cardiac arrhythmia diagnosis in the outpatient clinical setting. It includes a lightweight, cable-free adhesive Sensor that patients wear up to 7 days during their normal activities. Acquired Sensor ECG data is transferred via a user-supplied Windows® PC that connects to the Cardea SOLO Smart Cable. Cardea SOLO Software proprietary data processing enables onsite comprehensive arrhythmia data analysis, report generation, editing, storage and file sharing.

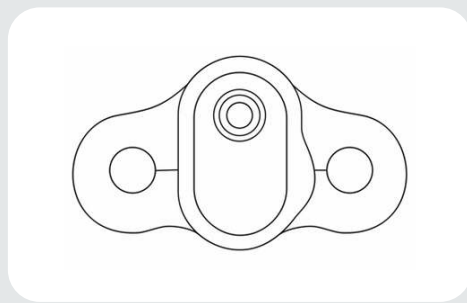

Sensor

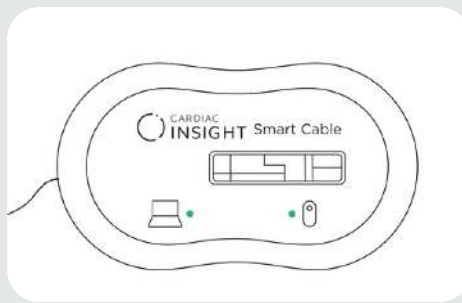

Smart Cable



Software Generated Report

Supported Operating System and Output Specifications

Operating System	Windows® 10 and higher
Software Installation	USB flash drive or electronic distribution.
Average File Size per Patient Test*	450 MB
Software Data Output Types	PDF, DAT, Excel files; one Windows folder per patient test. PHI/Patient demographics entered into PC and stored separately as an Excel .csv file. <i>No Protected Health Information (PHI) is acquired or stored in the single-use Cardea SOLO Sensor.</i>

*Cardea SOLO Software checks the available disk space for saving patient ECGs and associated information at start-up. If the available disk space is less than 10 gigabyte (GB) a message will be displayed.

Cardea SOLO Software

USER-SUPPLIED HARDWARE REQUIREMENTS

Computing	Windows® compatible personal computer
Hard disk	10 GB of free disk space or greater
Processor type and speed	Dual Core CPU @ 2.5 GHz or greater, 64-bit (x64) processor or equivalent
Display	1300 × 768 or higher resolution
System memory (RAM)	Minimum 8 GB
Pointing device	Windows® compatible pointing device
Keyboard	Windows® compatible keyboard
Ports	One available High Speed USB 2.0 port (minimum) for Smart Cable connection

USER-SUPPLIED SOFTWARE REQUIREMENTS

Network functionality/access	Per user system architecture. Patient data store on network drive.
EMR connectivity	EMR can access patient reports if stored on a shared network drive.
Remote access requirements services offered by the manufacturer	None
Security/password protection	Administrator and User-level password options.
EULA (End User Software License Agreement)	Displayed during software installation; EULA documentation available via "Help About."
Product support	Call (866) 554-3751, Option 2 from 6:00 AM to 5:00 PM (Pacific Time), Monday through Friday. Email: support@cardiacinsightinc.com .

Microsoft and Windows are trademarks of Microsoft Corporation.

Cardea SOLO® is a registered trademark of Cardiac Insight, Inc.

All product names, logos, and brands are property of their respective owners. All company, product and service names used in this document are for identification purposes only. Use of these names, logos, and brands does not imply endorsement. Refer to Cardea SOLO Operator's Manual and other product documentation for additional information and specifications.

Cardiac Insight reserves the right to make changes in specifications and features shown herein at any time and without notice or obligation.